

Combined push-button programme shift & volume control

Your hearing instrument has an easy to locate push-button for volume adjustment and switching between different programmes.

Programme shift:

Press button for 2 seconds and let go - the hearing aid will then beep twice indicating you are in program 2.


Prog.	No. of beeps	*Prog. assigned by audiologist
1	One	
2	Two	
3	Three	
4	Four	

* To be filled out by your audiologist


Volume adjustment:

Press briefly up or down to increase or decrease volume respectively.

Your hearing instrument will beep once to let you know volume has changed.

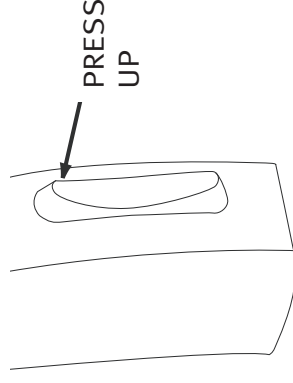

To move up

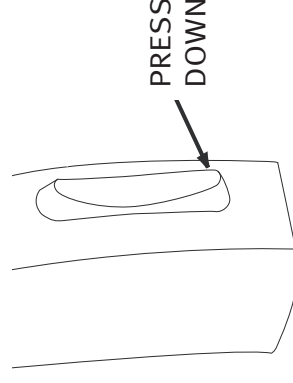

To move down

Additional options set by audiologist

Option	Volume Control	No. of short beeps	*Prog. assigned
1	Reaches preferred level	2	
2	Reaches Min & Max levels	3	

* Audiologist - please tick if assigned